

Doctor \_\_\_\_\_

Patient \_\_\_\_\_ Date \_\_\_\_\_

Street \_\_\_\_\_

Phone \_\_\_\_\_ Due Date \_\_\_\_\_

City \_\_\_\_\_ State \_\_\_\_\_ Zip \_\_\_\_\_

Return  UPS Gnd.  2nd day  Overnight  Other

**Please circle the appropriate items**

**Items shipped with this case**

Counter Model	Master Impression	Shade Tab	Bite Registration	Metal Articulator	Study Model
---------------	-------------------	-----------	-------------------	-------------------	-------------

**Teeth to Be Restored**

1 2 3 4 5 6 7 8 9 10 11 12 13 14 15 16  
 32 31 30 29 28 27 26 25 24 23 22 21 20 19 18 17




<b><u>Restoration Type:</u></b>	<i>Full cast</i>	<i>PFM</i>	<i>Captex</i>	<i>Empress</i>	<i>IPS Eris</i>
<i>Sculpture</i>	<i>Feldspathic Veneer</i>	<i>Fortress</i>	<i>Zirconia with pressed glass</i>	<i>Zirconia with Feldspathic</i>	<i>Milled Inceram</i>

<b><u>Please Return:</u></b>	Single units	Splint	Completed	Frame only	Bis Bake
------------------------------	--------------	--------	-----------	------------	----------

<b><u>Materials:</u></b>	50% gold	88% gold	Synspar Porcelain	IPS dSign Porcelain
--------------------------	----------	----------	-------------------	---------------------

<b><u>Templates:</u></b>	Waxup	Template and prep guide	Prepare teeth on duplicate model	Imaging date _____
--------------------------	-------	-------------------------	----------------------------------	--------------------

<b><u>Texture:</u></b>	Heavily defined	Moderate	Light	Smooth
------------------------	-----------------	----------	-------	--------

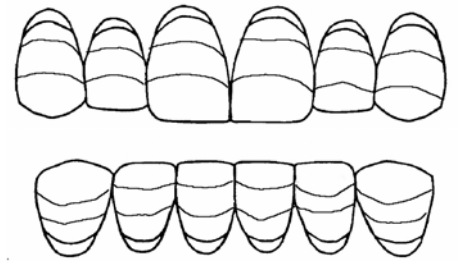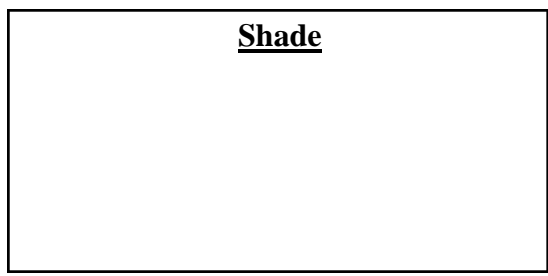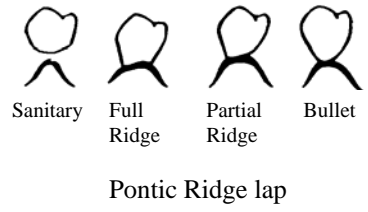

Additional Instructions: ▼

---

---

---

---

---

---

---

---

---

---

I want a photo preview of my case Y N Dr. Signature \_\_\_\_\_ Date \_\_\_\_\_

# Turn Your Imaging Case into a Wax Replica!

## The Template Technique

The Template Technique is a case-delivery protocol which **insures a predictable result before starting a cosmetic case**. Our dental technicians create a three-dimensional resin replica using cosmetic imaging as a guide. Then, we make a preparation guide and two-part temporary template from the replica. The preparation guide insures the accurate tooth reduction necessary to achieve an excellent result. The temporary template allows the dental team to fabricate temporaries in a fraction of the time it takes using other methods. In addition, **the temporaries look just like the case shown in the cosmetic imaging**. This new generation of temporaries offers patients an accurate preview of their final case. Once the temporary restorations are accepted by patients, our laboratory can fabricate permanent restorations – always according to plan. There are **no surprises!**

## Suggested steps for the Template Technique

1. **Verify with the patient that the cosmetic imaging pictures give an accurate representation of the result you are planning to deliver!** Never assume that patients “understand” anything they do not see. For example, never present Smile-Vision imaging which shows ten teeth treated if you plan to treat only four or six teeth.

2. **Send study models of your case to us at Smile-Vision.** By taking measurements from the cosmetic imaging pictures we create a 3-D resin replica of the planned result, a tooth preparation guide and a template for temporaries. (We will ask you to **Preview** images of the waxup on your imaging site before shipping).

3. **During reduction of the teeth**, insert a periodontal probe into the holes of the prep guide to verify correct preparation.

4. **After preparation**, follow these steps to use the two-part temporary template:

1. Insert the soft inner lining into the hard outer shell and try it over the prepared teeth to verify fit.
2. Coat the teeth with dentinal bonding agent or desensitizer.
3. Inject Bis Acrylic temporary resin into the template (Protemp 3 Garant by 3M—ESPE is our favorite).
4. Seat the template over the prepared teeth and apply considerable force while the resin cures on the teeth.
5. Remove the hard outer shell. Then peel the inner shell off the resin. The temporary restoration will have a smooth polished surface. You may reuse the template.
6. Trim the excess resin with a safe-end finishing bur. Adjust occlusion when necessary.
7. Verify patient satisfaction with the esthetic result 24 to 48 hours after placing the temporaries.

5. **For the final restorations** take an impression of the temporaries and send them to Smile-Vision for fabrication.



Resin Replica



Tooth preparation guide



Temporary template



Finishing temps

**Ask us to prep the model for you using the prep guide**

## Smile-Vision

143 California St. Newton, MA 02458 Telephone 800-634-3480 X 15  
Fax 617-924-8163 www.smilevision.net, e-mail info@smilevision.net